

A150 Benchtop Conductivity Meter

Instruction Manual

Introduction

Thank you for selecting the A150 benchtop conductivity meter. This manual provides a step-by-step guide to help you operate the meter, please carefully read the following instructions before use.

Unpacking

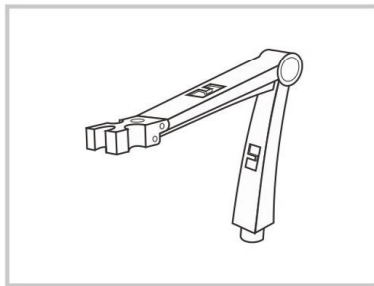
Before unpacking, ensure that the current work environment meets following conditions.

- Relative humidity is less than 80%.
- Ambient temperature is greater than 0°C and less than 60°C.
- No potential electromagnetic interference.

The following list describes the standard components of the meter. After the unpacking, please check all components are complete. If any are damaged or missing, please contact nearest distributor.




A150 Conductivity Meter




Electrode Holder



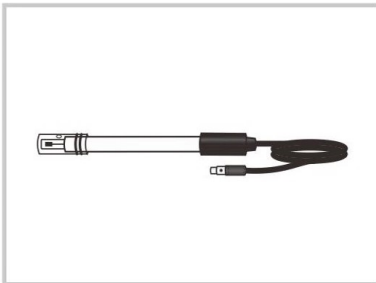
Conductivity Standard Solutions



CON-1 Conductivity Electrode

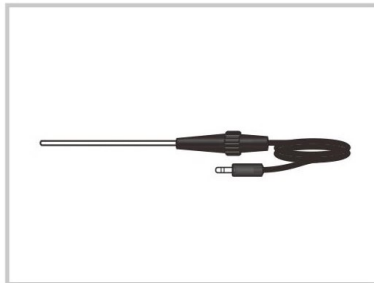



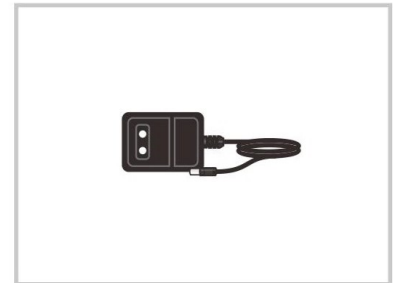





TP-10K Temperature Probe




DC12V Power Adapter


Keypad

The A150 meter has a succinct membrane keypad, names and symbols describe the each function key controls.

KEY	FUNCTION
	<ul style="list-style-type: none"> Switches the meter ON/OFF. Exits the calibration or setting and returns to measurement.
	<ul style="list-style-type: none"> Selects the measurement mode. Sets the temperature of sample (Press and hold the key for 3 seconds).
	<ul style="list-style-type: none"> Starts calibration. Enters the setup menu (Press and hold the key for 3 seconds).
	<ul style="list-style-type: none"> Locks the measured value. Resume measuring.
	<ul style="list-style-type: none"> Sends data to a printer or computer.
	<ul style="list-style-type: none"> Stores current reading to memory. Increase value or scroll up through the menu item.
	<ul style="list-style-type: none"> Views the calibration report or data logs. Decrease value or scroll down through the menu item.
	<ul style="list-style-type: none"> Confirms the calibration, settings or displayed options.

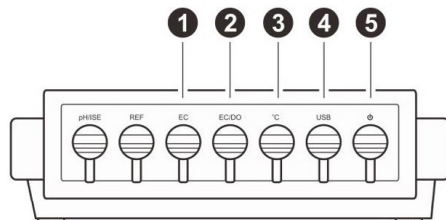
Connectors



NO.	CONNECTOR	DESCRIPTION
1	EC	Used for connecting the 4-pole conductivity electrode
2	EC/DO	Used for connecting the 2-pole conductivity electrode
3	°C	Used for connecting the temperature probe
4	USB	Used for connecting the computer or printer
5		Used for connecting the power adapter

Installing the Electrode Holder


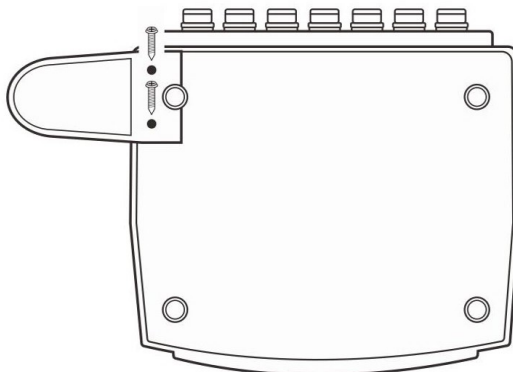
Take out the electrode holder from the packaging. Turn the meter over. Align the base plate of the electrode holder with the circular holes on the meter. Moderately tighten two screws.



Adjustment of electrode arm


After installation, if the electrode arm automatically rises or falls, you need to adjust the screw until arm locate at any position.

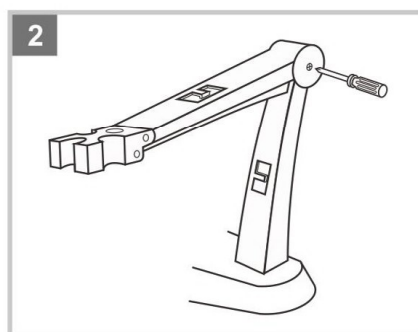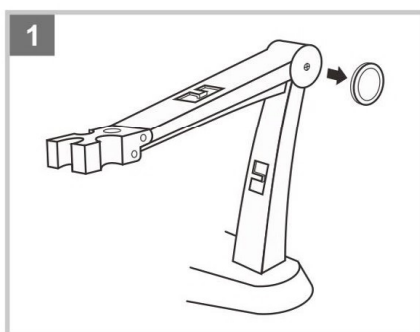1. Remove the plastic cover from the electrode arm.
2. Use the screwdriver to tighten the screw moderately.
3. Insert the plastic cover to previous position. Installation is completed.




Connecting the Electrode

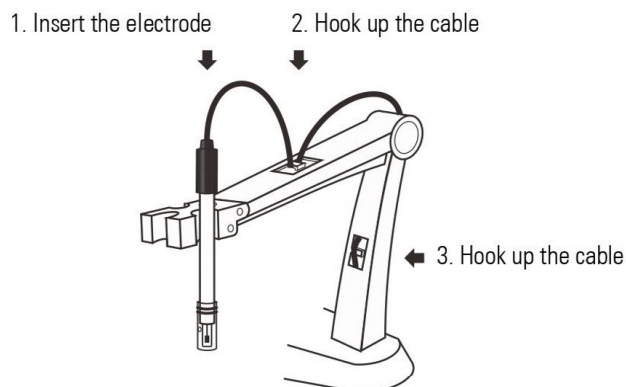1. Take out the conductivity electrode from the packaging. Follow the steps below to place the electrode into left or right side of the electrode arm.




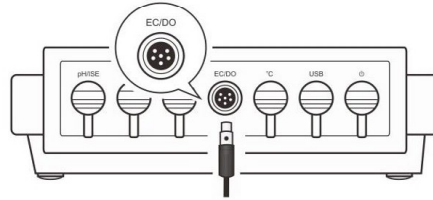2. Insert the 6-pin connector into the connector socket labeled EC/DO. After the connection is completed, DO NOT pull on the cable. Always make sure that the connector is clean and dry.




Connecting the Temperature Probe

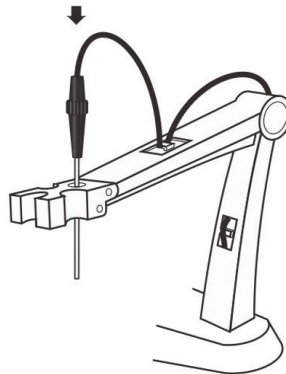1. Place the temperature probe into the circular hole of the electrode arm.




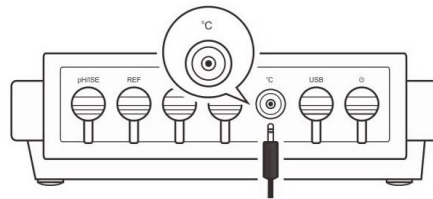2. Insert the phone plug to the connector socket labeled °C. Ensure the connector is fully seated.





Connecting the Power Adapter

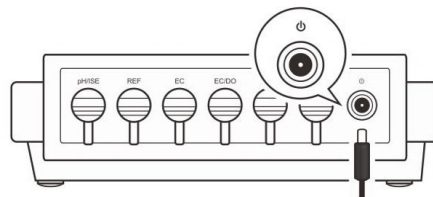1. Before plugging in the power adapter, ensure that its voltage matches the local main voltage.
2. Insert the connector to the power socket. The meter is now ready for use.



Switching the Meter On and Off

- Press and hold the  key to switch on the meter, the display shows the measured values.
- Press and hold the  key for 3 seconds, the meter will switch off.

 To enable the Auto-Power Off feature, please refer to chapter SETUP MENU.


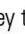

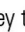

Setup Menu

The A150 conductivity meter contains an integrated setup menu that is used to customize the displayed option to meet measurement requirements. The following table describes the functions of the menu items.

MENU	OPTIONS	DESCRIPTION	DEFAULT
Sample ID	0000 to 9999	Set the sample ID to associate readings with the data log.	0000
Cell Constant	2-pole electrode (K=0.1)	Set the cell constant to match connected electrode.	2-pole electrode (K=1)
	2-pole electrode (K=1)		
	2-pole electrode (K=10)		
	4-pole electrode		
Calibration Points	1 to 3 points	Set the number of calibration points.	3 points
Temperature Coefficient	Linear (Range: 0.0 to 10.0%/°C)	Set the temperature compensation type and coefficient.	Linear (2.1%/°C)
	Non-linear		
Pure Water Coefficient	Enable	When the option is enabled, pure water coefficient will be applied automatically for ultra-pure water measurements.	Disable
	Disable		
Reference Temperature	20°C	Set the normalization temperature for conductivity measurement and calibration.	25°C
	25°C		
TDS Factor	Range: 0.01 to 1.00	Set the default TDS conversion factor.	0.50
Alarm Limits	Enable	Set the high and low limit values to activate alarm (Range: 0 to 999).	Disable
	Disable		
Calibration Due	Enable	Set the calibration interval to activate alarm (1 to 31 days).	Disable
	Disable		
Temperature Unit	°C	Set the default temperature unit.	°C
	°F		
Stability Criteria	Standard	Set when a measurement is recognized as stable.	Standard
	High-accuracy		
Auto-Read	Enable	When the option is enabled, the meter will automatically sense a stable reading and lock the measurements.	Disable
	Disable		
Auto-Power Off	Enable	When the option is enabled, the meter will automatically switch off if no key is pressed within 3 hours.	Disable
	Disable		
Date and Time	Year-month-day, hour-minutes	Set the current date and time.	

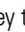

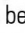
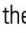
Interval Readings	Off	When the option is enabled, the meter will automatically send the measured data to the computer or printer.	Off
	10 seconds		
	30 seconds		
	60 seconds		
	10 minutes		
	30 minutes		
Password	Enable	Set the password protection for calibration and settings.	Disable
	Disable		
Brightness	Low, Mid, High	Set the brightness level of the backlight.	Mid
Clear Stored Data	Enable	Delete all stored readings in the memory.	Disable
	Disable		
Factory Reset	Enable	Reset the meter to factory default settings.	Disable
	Disable		






Setting the default option

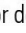1. In the measurement mode, press and hold the  key for 3 seconds to enter the setup menu.
2. Press the  or  key to select the menu item.
3. Press the **Enter** key, the cursor changes to highlight.
4. Press the  or  key to select the desired option.
5. Press the **Enter** key to confirm, the meter returns to the measurement mode. Setting is completed.


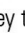



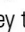
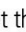




Setting the default parameter


The meter provides two methods for parameter settings.

- Press the  or  key to modify the value, press the **Enter** key to confirm.
- If the cursor appears below the first digit, press the  or  key to set the value, press the **Enter** key to confirm and move to the next digit. Repeat the steps above until the meter returns to the measurement mode. Setting is completed.

-  During the setting process, press the  or  key once, the setting value will increase or decrease gradually. Press and hold the  or  key, the setting value will increase or decrease quickly.

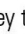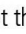Setup example - alarm limits

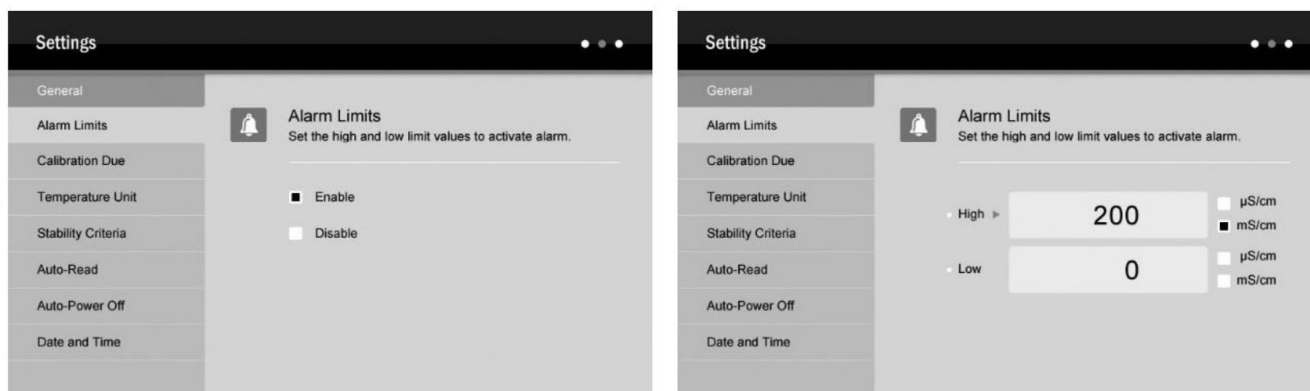
1. In the measurement mode, press and hold the  key for 3 seconds to enter the setup menu.
2. Press the  or  key to select the "Alarm Limits".
3. Press the **Enter** key, the cursor changes to highlight.
4. Press the  or  key to select the measurement unit, press the **Enter** key to confirm.
5. Press the  or  key to set the high alarm value, press the **Enter** key to confirm.
6. Press the  or  key to select the measurement unit, press the **Enter** key to confirm.
7. Press the  or  key to set the low alarm value, press the **Enter** key to return to the measurement mode.



Temperature Coefficient


The A150 conductivity meter contains two temperature correction methods. The linear temperature correction is appropriate for most samples. If the current sample is belong to the natural water (e.g., natural ground, well, or surface waters), using the non-linear correction is necessary. Note: non-linear correction can only be performed at temperature range from 0°C to 36°C. If the temperature value is out of above range, the meter will show a warning.

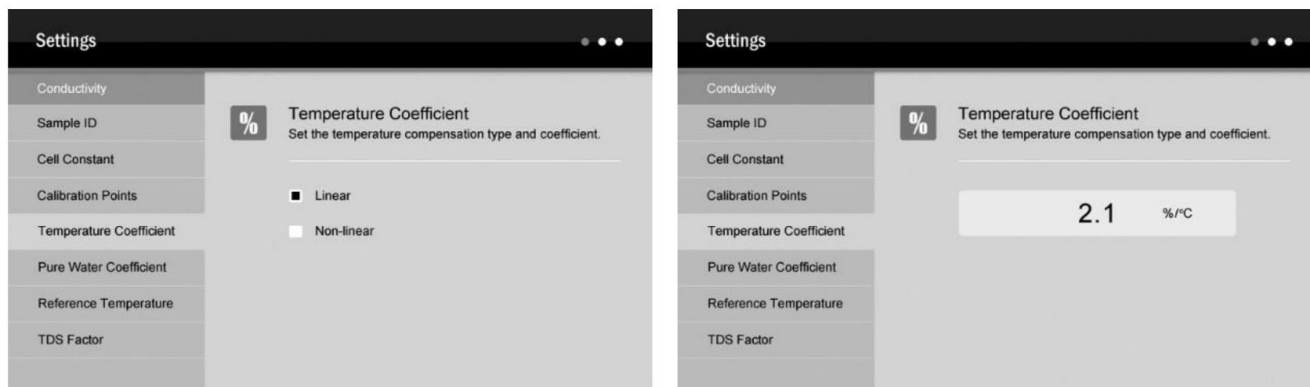❗ To calculate the temperature coefficient, please refer to addendum “How to calculate the temperature coefficient”.



Pure Water Coefficient

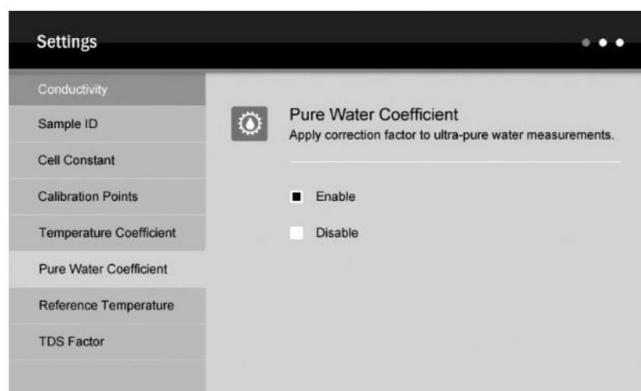The pure water coefficient is used to correct the samples with a conductivity of less than 5 μ S/cm. If enabled, the meter will be automatically calculated and applied coefficient for ultra-pure water measurements.




Stability Criteria

The Stability setting allows the user to set when a measurement is recognized as stable by the meter. When the Standard option is enabled, the Stable icon will quickly appear on the display. When the High-accuracy option is enabled, the icon will take longer to appear, but guarantees high accuracy of the measurement.




Auto-Read

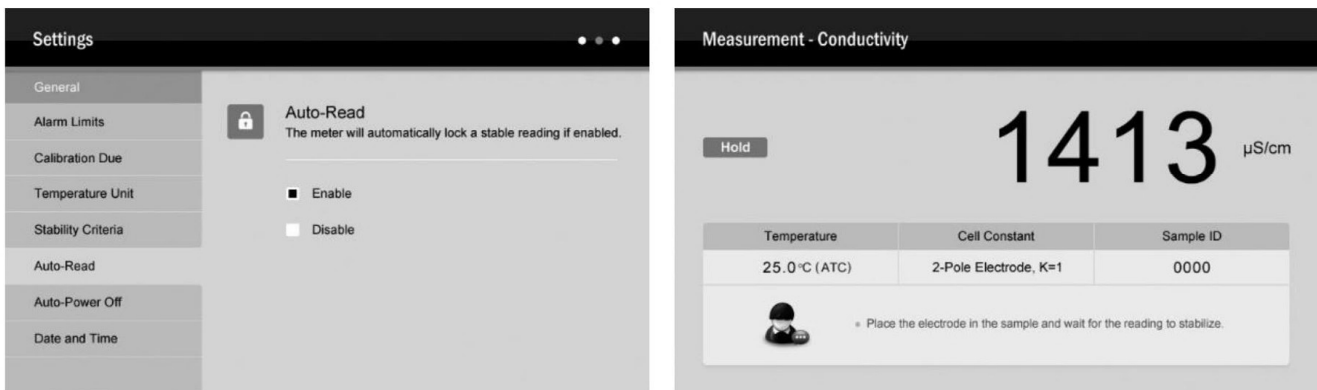The Auto-Read feature is used to lock a measurement endpoint. If enabled, the meter will automatically sense a stable reading and lock the measurements. The HOLD icon appears on the display. Press the **Meas** key, the meter resumes measuring.




Interval Readings

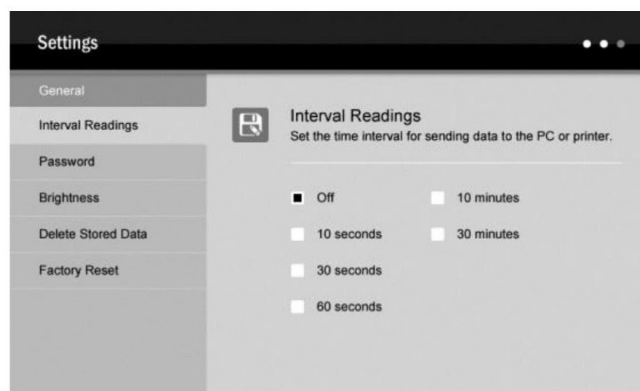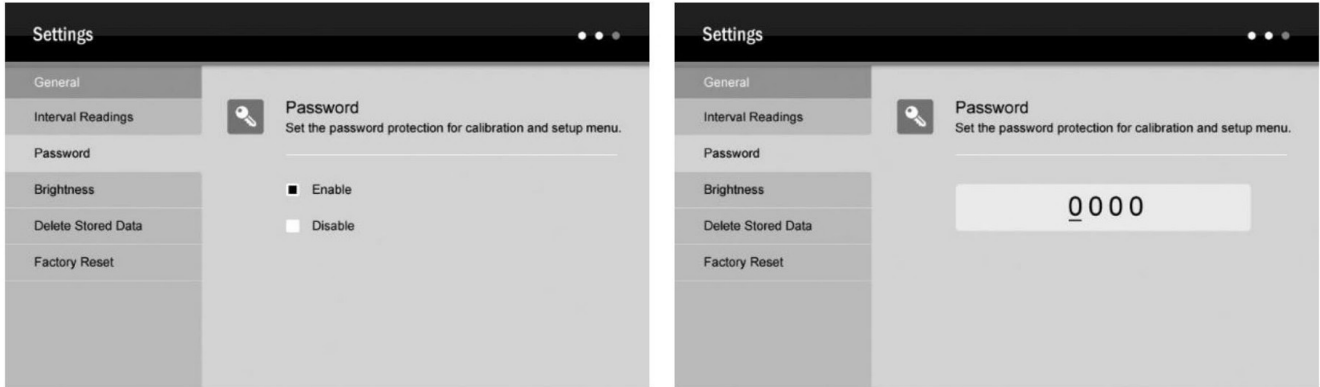The Interval Readings is capable of recording the measurements at the predefined time intervals. If enabled, the meter will continue to send measured data to the printer or computer until the measurement mode is exited. You are able to use the DAS software for receiving the data or viewing the real-time graph. For more details, please refer to chapter COMMUNICATION.






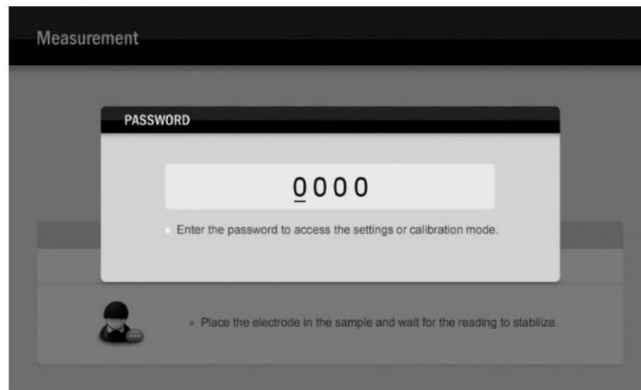
Password

The password protection is used to prevent the unauthorized calibration and settings. If enabled, the user must enter the 4-digit password to access the calibration or setup menu. If the setting value is 0000, the password protection will be invalid.



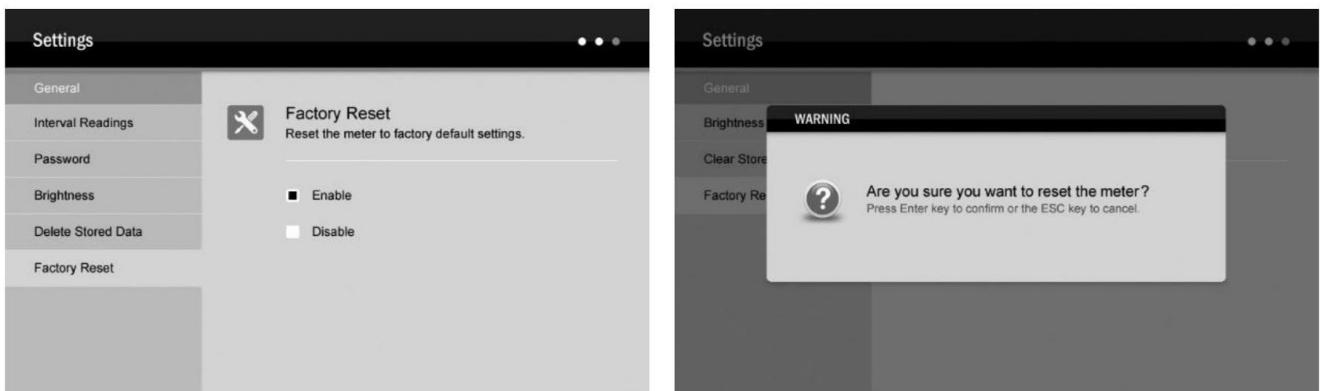
Unlock or Reset the password

Press and hold the  key in the measurement mode, the password input window immediately shows on the display and wait for entering the correct digits. Press the  or  key to input the password, press the **Enter** key to confirm. Once you have successfully entered the setup menu, selecting the "Disable". The password will be removed.



Factory Reset

The Factory Reset will restore the meter back to factory default settings. If enabled, all of the calibration data and selected options/parameters will be lost or reset, the meter must be recalibrated. During the setting process, when the display shows "Are you sure you want to reset the meter?", press the **Enter** key, the meter will immediately restore the factory settings, press the **ESC** key to cancel.



Prior to Use

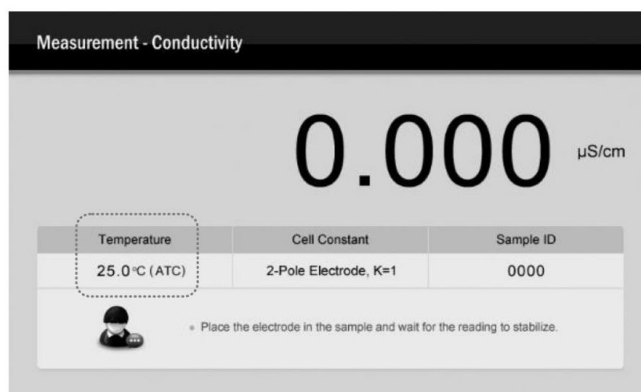
Soak the conductivity electrode for a few minutes in tap water to remove dirt and oil stains on the electrode.

Temperature Compensation

For better accuracy, we recommend the use of either a sensor with a built-in or a separate temperature probe for the calibration or measurement.

Automatic Temperature Compensation


Connect the temperature probe to the meter (Refer to page 4 "Connecting the Temperature Probe"). The ATC icon immediately appears on the display, the meter is now switched to the automatic temperature compensation mode.



Manual Temperature Compensation

If the meter does not detect a temperature probe, the MTC icon will show on the display indicating that the meter is switched to the manual temperature compensation mode. To set the temperature value, follow the steps below.

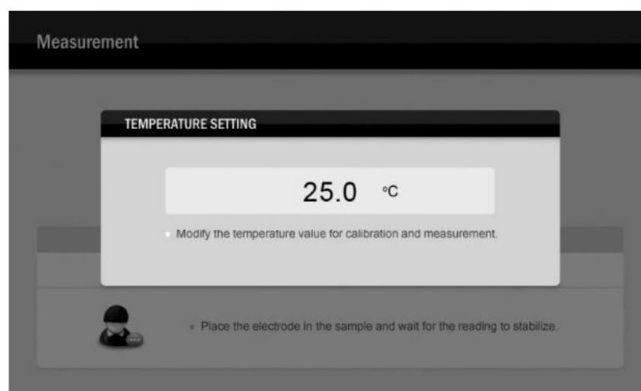1. Press and hold the °C key for 3 seconds to enter the temperature setting mode.
2. Press the ▲ or ▼ key to modify the temperature value.
3. Press the **Enter** key to confirm, the meter returns to the measurement mode. Setting is completed.



- ① Press the ▲ or ▼ key once, the setting value will increase or decrease by 0.1. Press and hold the ▲ or ▼ key, the setting value will increase or decrease by 1.

Selecting the Conductivity Electrode

The A150 conductivity meter is capable of using the 3 types of 2-pole conductivity electrodes. Prior to the calibration and measurement, you need to select the different sensors as per your sample concentration. The table below shows available electrodes and its effective measuring ranges.

CONDUCTIVITY ELECTRODE	MEASURING RANGES	CELL CONSTANT
CON-0.1	0.1~100 μ S/cm	K=0.1
CON-1	10 μ S/cm~10mS/cm	K=1
CON-10	100 μ S/cm~200mS/cm	K=10

If the 4-pole conductivity electrode is selected, its best measurement range is 100 μ S/cm to 200mS/cm.

Conductivity Calibration


The meter allows 1 to 3 points calibration in the conductivity mode. To ensure higher accuracy, we recommend that you perform 3 points calibration or select a standard value close to the sample value you are measuring. The meter will automatically detect these conductivity standard solutions and prompt the user to calibrate the meter. When the calibration is completed, all new calibration values will automatically override existing data.

The following table shows acceptable conductivity range of standard solution for each measuring range.


MEASURING RANGE	CALIBRATION SOLUTION RANGE	DEFAULT
0~20 μ S/cm	7~17 μ S/cm	10 μ S/cm
20~200 μ S/cm	70~170 μ S/cm	84 μ S/cm
200~2000 μ S/cm	700~1700 μ S/cm	1413 μ S/cm
2~20mS/cm	7~17mS/cm	12.88mS/cm
20~200mS/cm	70~170mS/cm	111.8mS/cm

Single point calibration

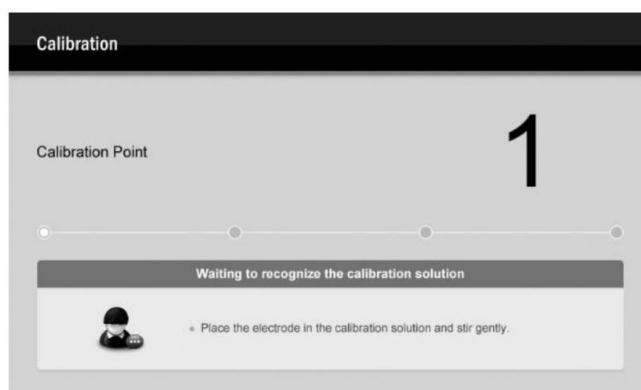
- 1.1 Ensure that the meter is in the conductivity measurement mode and you have selected 1 point calibration in the setup menu.
- 1.2 Rinse the conductivity electrode with distilled water, then rinse with a small amount of standard solution.
- 1.3 Press the **Cal** key, the display shows "Calibration Point 1" and waits for recognizing the standard solution.



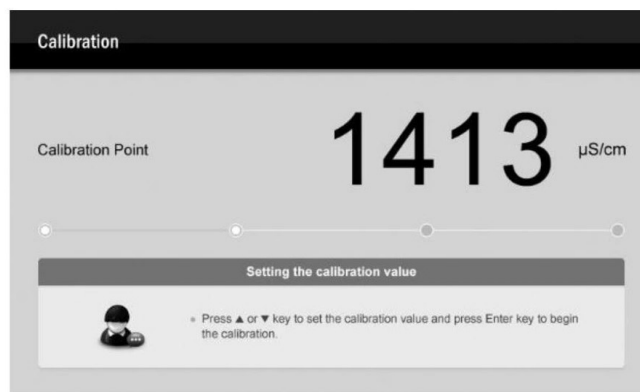
- 1.4 Place the conductivity electrode into the standard solution, the meter automatically shows current calibration standard (e.g., 1413 μ S/cm).



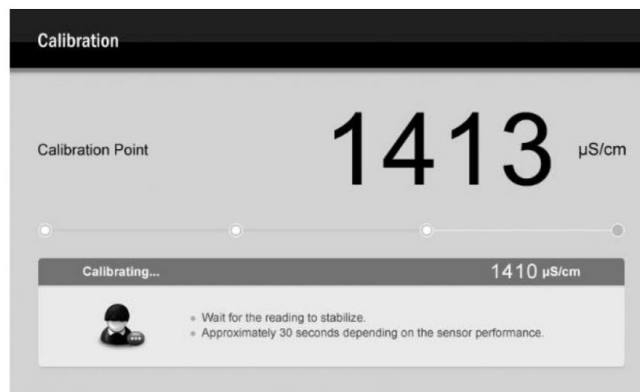- 1.5 If necessary, press the ▲ or ▼ key to modify the calibration value.
1.6 Press the **Enter** key, the Calibrating.. icon shows on the display.




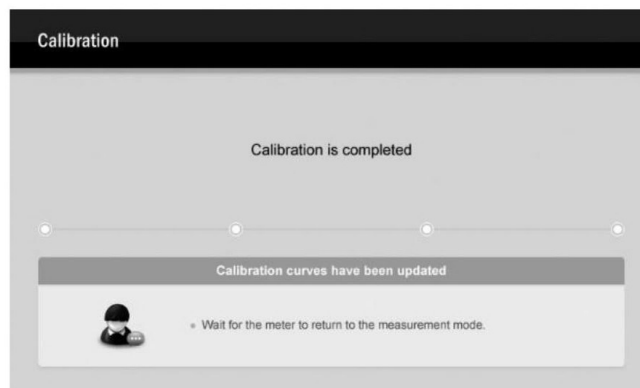- 1.7 Wait for the reading to stabilize, the meter automatically shows "Calibration is completed" and returns to the measurement mode.




Multi-point calibration

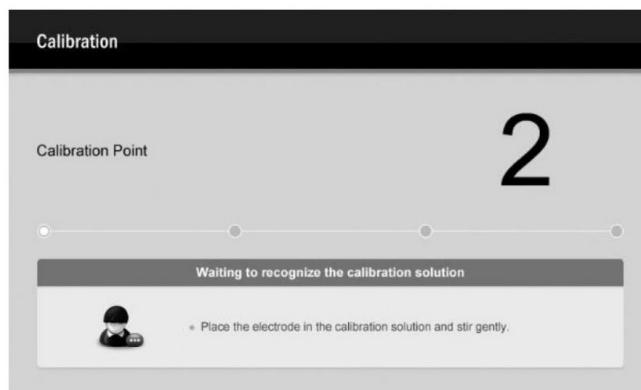
- 2.1 Ensure that you have selected 2 or 3 points calibration in the setup menu.
- 2.2 Repeat steps 1.2 to 1.6 above. When the first calibration point is completed, the display will show "Calibration Point 2". The meter prompts you to continue with second point calibration.



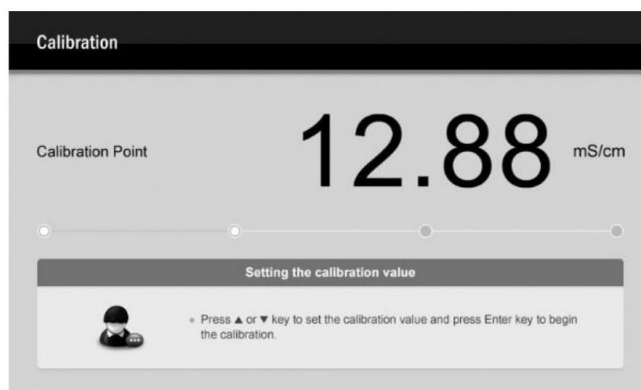
- 2.3 Rinse the conductivity electrode with distilled water. Place the electrode into the next standard solution, the meter automatically shows current calibration standard (e.g., 12.88mS/cm).



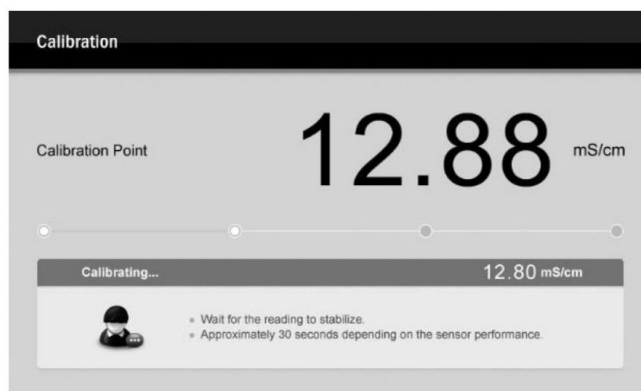- 2.4 If necessary, press the ▲ or ▼ key to modify the calibration value.
- 2.5 Press the **Enter** key, the meter begins the calibration.



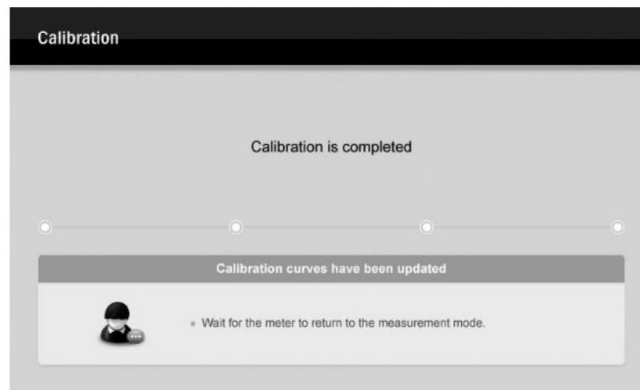- 2.6 Wait for the reading to stabilize, the display will show "Calibration Point 3". The meter prompts you to continue with third point calibration.
- 2.7 Repeat the steps 2.3 to 2.5 above until the meter returns to the measurement mode. Calibration is completed.




- i**
- Performing the conductivity calibration will simultaneously calibrate the corresponding TDS, salinity, resistivity and conductivity ash value.
 - If you want to exit the calibration, press the **ESC** key, the meter will immediately return to the measurement mode.

Temperature Calibration

During the measurement process, if the temperature reading displayed differs from that of an accurate thermometer, the meter needs to be calibrated.

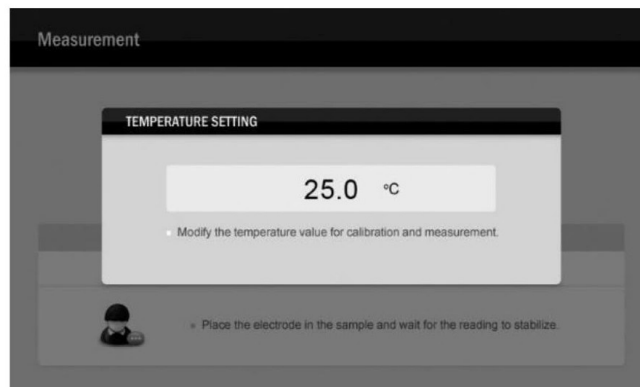
1. Connect the temperature probe to the meter and place into a solution with a known accurate temperature.
2. Press and hold the **°C** key for 3 seconds to enter the calibration mode.
3. Press the **▲** or **▼** key to set the temperature value.
4. Press the **Enter** key, the meter returns to the measurement mode. Calibrating is completed.



Calibration Report

The A150 conductivity meter provides detailed report for the conductivity calibration.

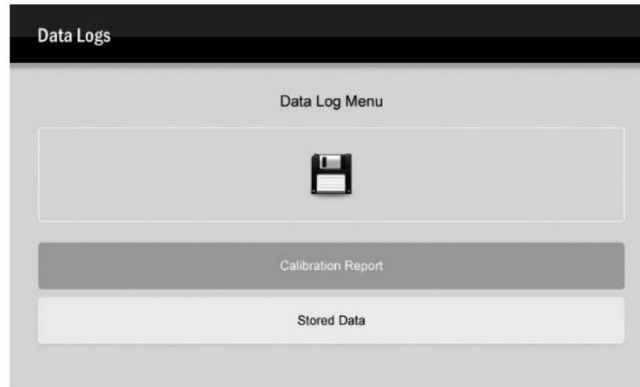1. Press the **MR** key in the measurement mode, the meter shows the data log options.




2. Press the **▲** or **▼** key to select the "Calibration Report".
3. Press the **Enter** key, the display shows the updated calibration information.
4. Press the **ESC** key, the meter returns to the measurement mode.

2-Pole Electrode			
Date:	2018 - 1 - 10	Reference Temp.:	25 °C
Time:	10 : 15 : 58	Calibration Due:	1 Day(s)
Calibration Points	Constant (K=1)	Constant (K=10)	Constant (K=0.1)
10.00 µS/cm	1.00	10.00	0.100
84.0 µS/cm	1.00	10.00	0.100
1413 µS/cm	0.98	10.00	0.100
12.88 mS/cm	0.98	10.00	0.100
111.8 mS/cm	1.00	10.00	0.100

Conductivity/TDS/Salinity/Resistivity Measurement

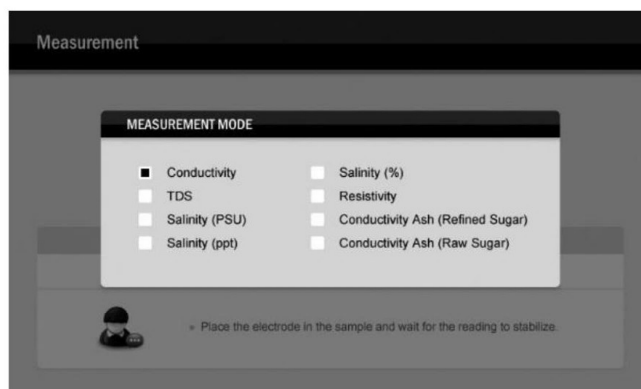
1. Press the **Mode** key in measurement mode and the **▲** or **▼** key to select desired option, press the **Enter** key to confirm.
2. Rinse the conductivity electrode thoroughly with distilled water.
3. Place the electrode into the sample solution, stir the electrode gently.
4. Record the measured value when the reading is stable.



Conductivity Ash Measurement

The A150 conductivity meter contains two conductivity ash measurement modes - Refined Sugar (ICUMSA GS2/3-17 standard) and Raw Sugar (ICUMSA GS1/3/4/7/8-13 standard).

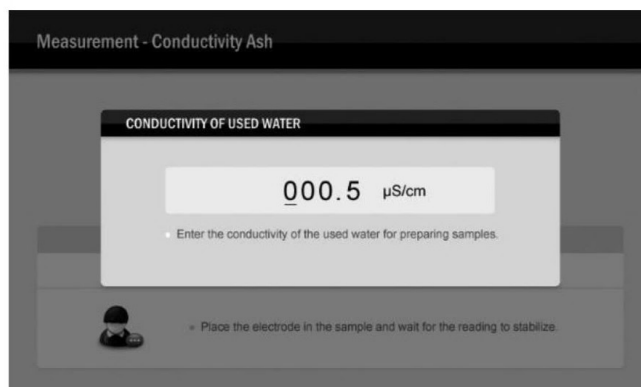1. Prepare the sugar sample according to the selected ICUMSA method.
2. Press the **Mode** key in measurement mode and the **▲** or **▼** key to select desired option, press the **Enter** key to confirm.
3. The display automatically shows a input window and waits for entering the conductivity of the used water for preparing sugar solutions (range: 0.0 to 100.0 $\mu\text{S}/\text{cm}$).



4. Press the **▲** or **▼** key to set the value, press the **Enter** key to confirm and move to the next digit. When the setting is completed, the meter will automatically enter the measurement mode.
5. Rinse the conductivity electrode thoroughly with distilled water.
6. Place the electrode into the sample solution, stir the electrode gently.
7. Record the measured value when the reading is stable.

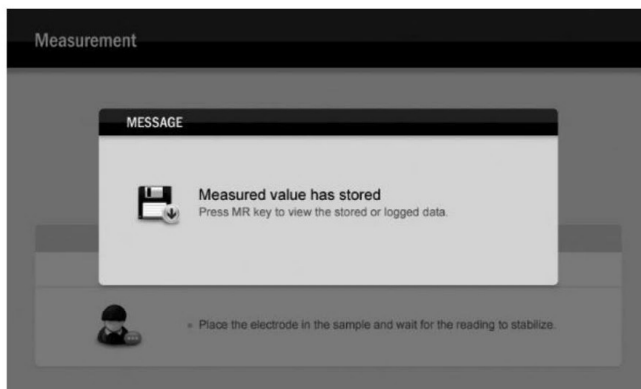
ⓘ NOTE: conductivity ash measurement can only be performed at temperature range from 15°C to 25°C. If the temperature value is out of this range, the meter will show a warning.

Storing and Recalling Data from Memory

The A150 conductivity meter are capable of storing and recalling up to 1000 data sets.


Storing a measurement result

During the measurement process, press the **MI** key to store the measured value, the meter will show a reminder as follow.



Recalling from memory

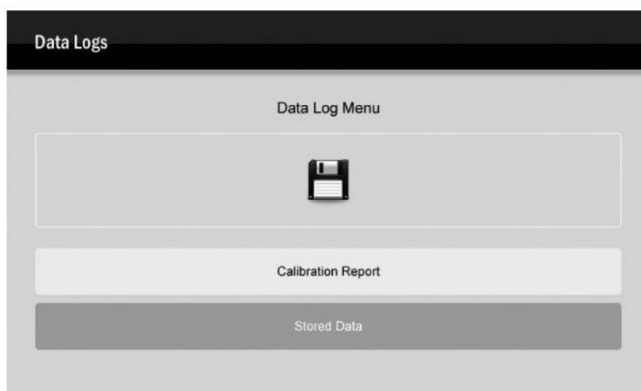- 1.1 Press the **MR** key in the measurement mode, the meter shows the data log options.
- 1.2 Press the **▲** or **▼** key to select the "Stored Data".




- 1.3 Press the **Enter** key, the display shows the data list.
- 1.4 Press the **ESC** key, the meter returns to the measurement mode.

Date and Time	Sample ID	Reading	Temperature
2018 - 1 - 10 10:12:30	0000	1413 μ S/cm	25°C
2018 - 1 - 10 10:12:55	0000	1308 μ S/cm	25°C
- - : :			
- - : :			
- - : :			
- - : :			
- - : :			



Clearing the memory

If the memory is full, the meter will automatically show a reminder and wait for user to delete all stored readings.

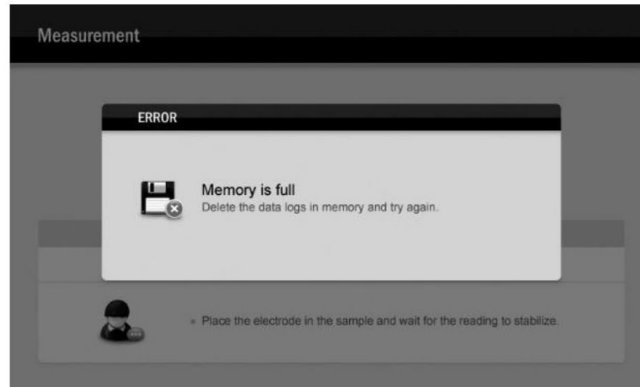WARNING: once the data are deleted that can not be recovered.




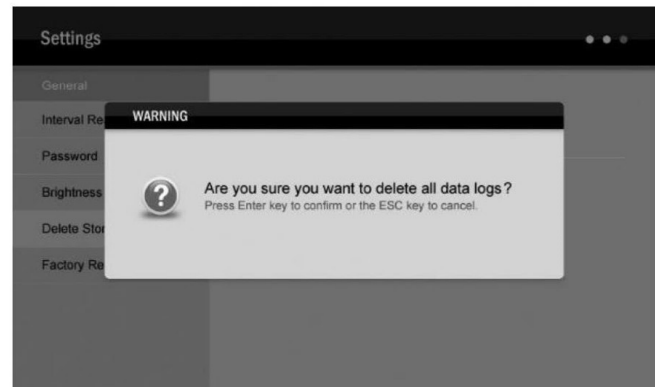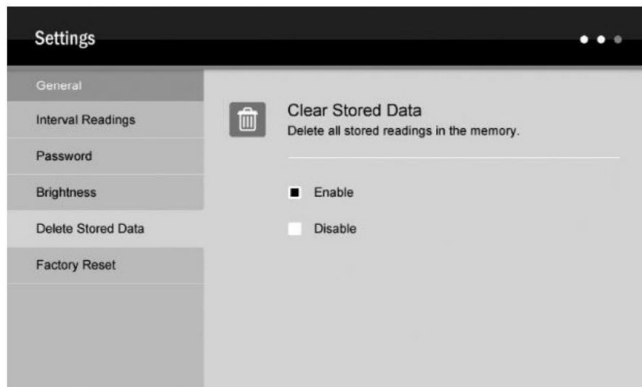- 2.1 In the measurement mode, press and hold the **⚙** key for 3 seconds to enter the setup menu.
- 2.2 Press the **▲** or **▼** key to select the "Clear Stored Data".
- 2.3 Press the **Enter** key, the cursor change to highlight.
- 2.4 Press the **▲** key to select the "Enable".
- 2.5 Press the **Enter** key, the meter shows a warning "Are you sure you want to delete all date logs ?"
- 2.6 Press the **Enter** key to confirm or the **ESC** key to cancel. The meter returns to the measurement mode.



Communication

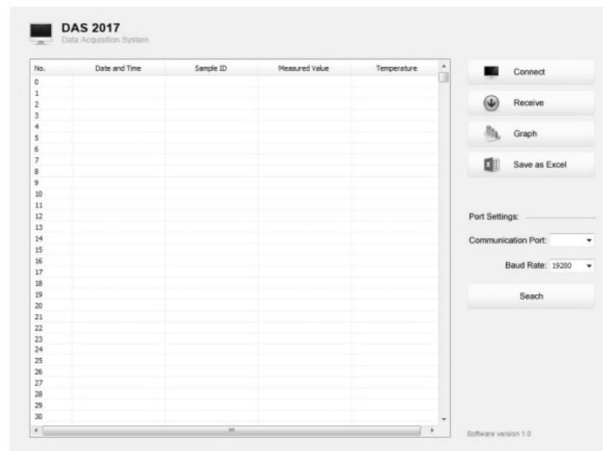Bante Instruments provides a Data Acquisition System that can be used to transfer data, receive the measuring values or import the data to Excel. You are able to download this software from our official website at www.bante-china.com. Before installation, ensure that Windows 7/8/10 operating system has been installed on your computer.



Receiving data

1. Connect the USB cable and data converter to the meter and computer.
2. Click the DAS_A_Series icon on computer, the system will automatically scan an available communication port and show the message box "Found a port on your computer".
3. Click the **OK** button, the application starts.
4. Click the **Connect** button, the screen shows "Port is connected" indicating that the communication between the meter and the computer has been established.
5. Click the **OK** button to confirm.
6. Click the **Receive** button, the stored data automatically transfer to computer.

Interval recording

This function is used to record the measuring value within the specify time period. The setting method refers to page 6 "Setting the Default Options". Note:

- The first data need 1 minute to be shown on screen.
- Do not press any key on meter during the Interval Recording mode, it will cause the communication interruption.

Graph mode

This function helps user to view variations of the measured value continuously. Click the **Graph** button, the screen immediately shows the curve graph. Click the **X** button to quit.

Create the excel file

When the data transfer is completed, click the **Save as Excel** button, the measured values in the data sheet will automatically convert to Excel file. **WARNING:** Once the software is closed, all received data will be lost and can not be recovered.

Conductivity Electrode Care and Maintenance

- DO NOT touch the measurement area of electrode and always make sure that is clean.
- If there is a build-up of solids inside the measurement area, these should be removed very carefully with a cotton bud soaked in solvent, taking care not to touch the metal parts of the inner cell. After the solids are removed, the electrode must be recalibrated.

Specifications

Conductivity	Model	A150
	Range	0.01~20.00, 200.0, 2000 μ S/cm, 20.00, 200.0mS/cm
	Accuracy	\pm 0.5% F.S
	Resolution	0.001, 0.01, 0.1, 1
	Calibration Points	1 to 3 points
	Calibration Solutions	10 μ S/cm, 84 μ S/cm, 1413 μ S/cm, 12.88mS/cm, 111.8mS/cm
TDS	Range	0.00mg/L~100.0g/L (Max. 200g/L)
	Accuracy	\pm 1% F.S
	TDS Factor	0.01~1.00 (Default 0.5)
Salinity	Range	0.00~80.00ppt, 0.00~42.00psu, 0.00~8.00%
	Accuracy	\pm 1% F.S
	Measurement Modes	Practical Salinity (psu), Natural Seawater (ppt) or %
Resistivity	Range	0.00~30.00M Ω
	Accuracy	\pm 1% F.S
Conductivity Ash	Range	0.00~100%
	Accuracy	\pm 1% F.S
	Measurement Modes	Refined Sugar or Raw Sugar
Temperature	Range	0~105°C, 32~221°F
	Accuracy	\pm 0.5°C
	Calibration Points	1 point
General	Temperature Compensation	0~100°C, 32~212°F, Manual or Automatic
	Temperature Coefficient	Linear, Non-linear and Pure water compensation
	Normalization Temperature	20°C or 25°C
	Cell Constant	2-pole electrodes (K=0.1, 1, 10) or 4-pole electrode
	Memory	Stores up to 1000 data sets
	Output	USB Communication Interface
	Power Requirements	DC12V/2A, using AC adapters, 220VAC/50Hz
	Dimensions	240 (L) \times 220 (W) \times 80 (H)mm
	Weight	1.7kg

Addendum 1: Preparation of Conductivity Standard Solutions

Place AR potassium chloride reagent in a 50ml beaker and dry in an oven for 3 to 5 hours at 105°C, then cool to room temperature in desiccator.

CALIBRATION SOLUTIONS	REAGENT
84µS/cm	Accurately weigh out 42.35mg of KCL and dissolve in 1 litre deionised water.
1413µS/cm	Accurately weigh out 745.9mg of KCL and dissolve in 1 litre deionised water.
12.88mS/cm	Accurately weigh out 7.45g of KCL and dissolve in 1 litre deionised water.
111.8mS/cm	Accurately weigh out 74.5g of KCL and dissolve in 1 litre deionised water.

Stir the solution until the reagent has thoroughly mixed. Preparation is completed.

Addendum 2: How to calculate the temperature coefficient

To determine the temperature coefficient of sample solution use the formula below:

$$T_c = \frac{C_{TB} - C_{TA}}{C_{TA}(T_B - 25) - C_{TB}(T_A - 25)} \times 100\%$$

Where:

T_c = Temperature coefficient

C_{TA} = Conductivity at Temperature A

C_{TB} = Conductivity at Temperature B

T_A = Temperature A

T_B = Temperature B

1. Press and hold the °C key for 3 seconds to enter the temperature setting.
2. Press the ▲ or ▼ key to set the temperature to 25°C.
3. Place the electrode into the sample solution A and record the temperature value T_A and conductivity value C_{TA} .
4. Condition the sample solution and electrode to a temperature that is about 5°C to 10°C different from T_A . Record the temperature value T_B and conductivity value C_{TB} .
5. Calculate the temperature coefficient according to the formula above.

Addendum 3: How to calculate the TDS conversion factor

To determine the TDS conversion factor use the formula below:

$$\text{Factor} = \frac{\text{Actual TDS}}{\text{Actual Conductivity @ 25°C}}$$

Where:

Actual TDS: value from the high purity water and precisely weighed NaCl or KCL reagent.

Actual Conductivity: the meter measured conductivity value.

For example: dissolve 64g of potassium chloride reagent in 1 litre distilled water. If its conductivity value is 100mS/cm, then TDS conversion factor is 0.64.

Hazardous Substance Statement

Instruments is committed to the reduction and eventual elimination of all hazardous substances in both the manufacturing process and finished products we supply. We have an active manufacturing and procurement program to minimize and eliminate the use of harmful heavy metals such as cadmium, lead, mercury and the like. New technologies and design parameters are also promoting these efforts and we expect to have little or no such materials in our product in the coming years. We welcome our customer suggestions on how to speed up these efforts.



Warranty

The warranty period for meter is one year from the date of shipment. Above warranty does not cover the sensor and calibration solutions. Out of warranty products will be repaired on a charged basis. The warranty on your meter shall not apply to defects resulting from:

- Improper or inadequate maintenance by customer.
- Unauthorized modification or misuse.
- Operation outside of the environment specifications of the products.

For more information, please contact the nearest authorized distributor.