


KD-2258 Rotary Microtome

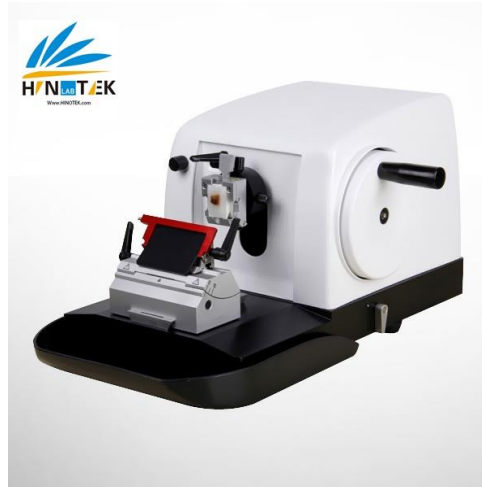

Features

- Ergonomic design, greatly reducing the discomfort associated with long hours of operation
- Ideal sectioning for using either a standard blade or disposable blade
- Cassette holder can be locked at any position along each of the three axes, enabling easy adjustment of specimen slicing angle
- Right operation hand-wheel can be locked at any position, providing maximum safety and convenience
- Specimen retraction and replaceable cassette holder system
- Blade holder base can be moved forward-or-backward and left-or-right flexibly to eliminating any contact of the fingers with the blade
- Convenient N-type blade holder and E-type holder for disposable blades
- The red bar on the blade holder covers the whole length of blade to protect the user and the push bar enables easy changing of the blade
- Easy-to-clean well-sealed housing does not need to be removed for cleaning
- Large-volume waste tray is easy to remove and clean

Technical Specifications

- Section Thickness Setting Range: 0 – 60µm
- Section thickness setting: 0 – 60µm
- 0 – 2µm increment 0.5µm
- 10 – 20µm increment 2µm
- 20 – 60µm increment 5µm

Precision Error	±1%
Minimum Setting of Sectioning Thickness	0.5µm
Maximum Specimen Size	60×50 mm
Total Horizontal Specimen Feed	24 mm
Vertical Specimen Stroke	52 mm
Specimen Retraction	0 – 28 mm
Trimming Thickness	1 – 60µm
Dimensions	565×440×285mm (W×D×H)
Net weight	27kg