

DLAB LED Digital Rotary Evaporator is a basic chemical laboratory instrument suitable for evaporation, distillation and separation of chemicals. It can be used as a multifunctional system in combination with vacuum pump, vacuum controller and circulating refrigerators to provide excellent distillation solutions in many laboratory applications.

Features

- 5L heating bath with a wide temperature range (room temp. to 180°C). Independent temperature control makes it can be used separately. Water/oil heating mode can be changed only through a switch
- Manual lift combined with auxiliary lift for precise positioning of the glassware
- PID control ensures high temperature accuracy at $\pm 1^\circ\text{C}$ (water)
- Overheating protection temperature at 220°C
- Boil-dry protection, automatically power off if heating without water/oil in the heating bath
- Speed range from 20 to 200rpm, timing interval operation in clockwise and anticlockwise directions for drying process
- Patented condenser (cooling surface 1200/1700cm² for selection), double helix condensing tube + center circular arc design to speed up the flow rate of liquid with excellent cooling effect
- Ejection mechanism ensures easy exchange of evaporating flask
- Adjustable immersion angle
- Patented double spring sealing ring which is made of PTFE provides an excellent sealing performance
- Remote function provides PC control and data transmission
- Compatible with the entire range of RE100-Pro glassware
- Optional glassware with explosion proof film available
- Overheating warning above 50°C



USB connector



Chemical resistance




Overheating protection



Compatible with all range of glassware

RE100-S LED Digital Rotary Evaporator

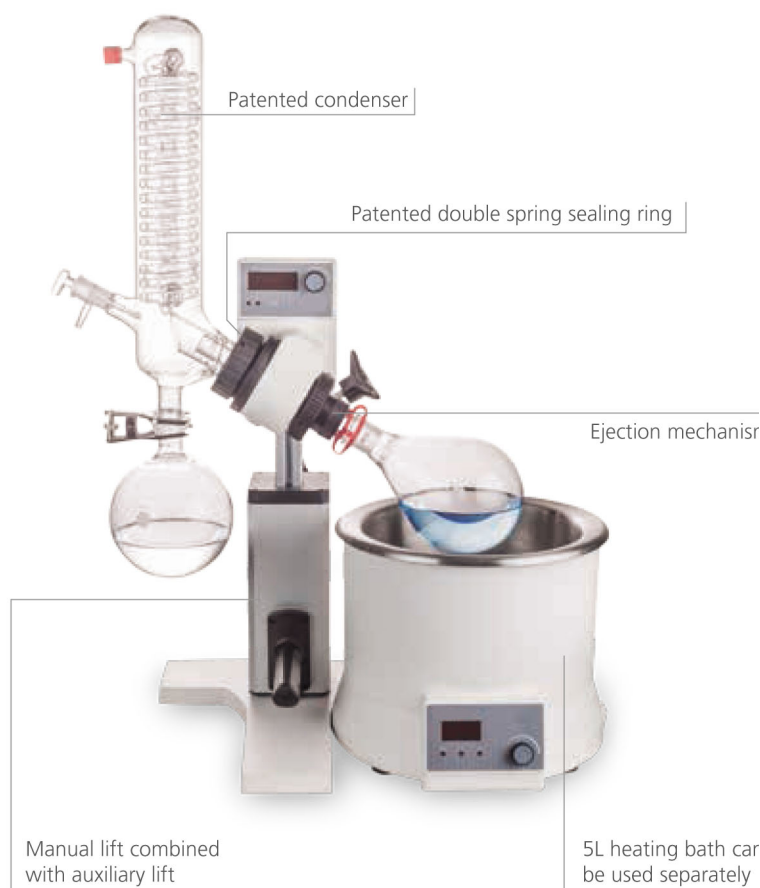

Specifications

RE100-S	
Motor Type	Brushless DC motor
Speed Range	20-200rpm
Display	LED (speed, temperature, time)
Clockwise and Anti-clockwise	Yes
Heating Temperature Range	Room temp. to 180°C
Control Accuracy	water: $\pm 1^\circ\text{C}$ oil: $\pm 3^\circ\text{C}$
Heating Power	1200 W
Stroke Displacement	Manual 110mm + auxiliary 100mm
Interval Time Setting Range	1~999s
Dimension [D×W×H]	Main Unit: 440×320×450mm Heating Bath: 300×300×240mm
Weight	Main Unit: 7kg Heating Bath: 3kg
Permissible Ambient Temperature	5~40°C
Permissible Relative Humidity	80%RH
Protection Class	IP20
USB Interface	Yes
Voltage/Frequency	100~120/200~240V 50/60 Hz
Power	1245 W

Condenser



New

Vertical

Cat. No.	Glassware type	Volume	Cooling surface
18900165	Evaporating flask(NS29/32) Receiving flask(KS35/20)	1000mL	1700cm ²
18900173	Evaporating flask(NS24/40) Receiving flask(KS35/20)	1000mL	1700cm ²
18900166	(COATED)Evaporating flask(NS29/32) Receiving flask (KS35/20)	1000mL	1700cm ²
18900174	(COATED)Evaporating flask(NS24/40) Receiving flask (KS35/20)	1000mL	1700cm ²
18202410	evaporating flask (NS 29/32) Receiving flask (KS 35/20)	1000mL	1200cm ²

Diagonal

Cat. No.	Glassware type	Volume	Cooling surface
18900163	Evaporating flask(NS29/32)Receiving flask(KS35/20)	1000mL	1700cm ²
18900171	Evaporating flask(NS24/40)Receiving flask(KS35/20)	1000mL	1700cm ²
18900164	(COATED)Evaporating flask(NS29/32) Receiving flask (KS 35/20)	1000mL	1700cm ²
18900172	(COATED)Evaporating flask(NS24/40) Receiving flask (KS35/20)	1000mL	1700cm ²

Heating Bath

Cat. No.	Description	Volume
18900201	Rotary Evaporator, 220V	5L
18900202	Rotary Evaporator, 110V	5L
6033100210	LED heating bath, 100V-120/200-240V,50/60HZ	5L

Heating Bath Lid

Cat. No.	Description
18900547	heating Bath lid for RE100-S, improve distillation efficiency & ensure personal safety from liquid splash

