


Product Name: Automatic Seed Counter

Model: SLY -E Plus

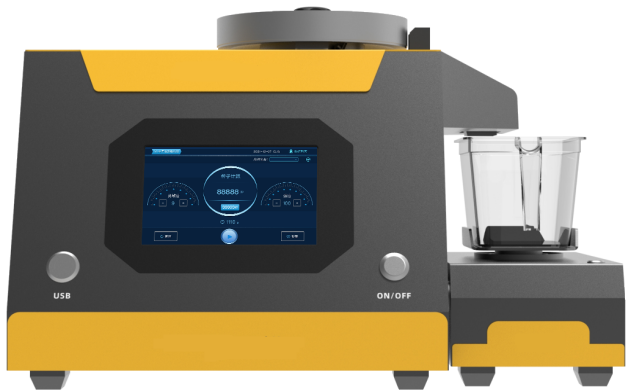

Characteristics

- Super-large touch screen design: 7-inch color touch screen, enhance the user's human-computer interaction experience
- Adjustable speed of counting disc: with stepless speed regulation function, the speed of counting grains can be adjusted
- Accurate counting: the discharge channel can be adjusted, only counting the particles of the selected magnitude, avoiding interference such as impurities
- High-precision weighing: built-in high-precision gravity sensor, can automatically weigh, automatically count grains, and automatically calculate thousand-grain weight
- Automatic deceleration, more accurate counting: decelerate and count after approaching the set number of seeds to improve the counting accuracy
- Automatic data recording: the measured data can be automatically recorded and saved, and the historical data can be viewed at any time
- Data export: support USB data export

Technical specification

| | |
|-----------------|--|
| Seed length | 0.7~23mm |
| Counting error | ±2‰ |
| Weighing error | ±0.2g |
| Vibration noise | ≤70dB |
| Counting speed | ≥1000 grains/3 minutes (Speed can be adjustable) |
| Weighing range | 1~1000.0g |

Related Products

Grain Moisture Meter (with volume-weight)/Grain Moisture Meter/Seed Blower/Seed Neatness Workbench